



**More personal. More plasma.**

**HAEMONETICS®**



More personal. More plasma.

## Donors come in every size.

Now, so do collection volumes, with the **Persona™ Plasma Collection Solution.**

**Persona™**, a proprietary integrated plasma collection solution built upon the **NexSys PCS® platform**, is the first and only donor-tailored solution clinically shown to yield +9% to 12%\* more plasma per donation on average to maximize both cost-efficient output and patient impact.

**+9% to 12%\*** AVERAGE INCREASE IN PLASMA VOLUME



\*Based on baseline device, software configuration and donor population

### THE PERSONA™ NOMOGRAM

The Persona Plasma Collection Solution uses the **Persona Nomogram**—a novel, personalized percent plasma nomogram that tailors donation targets to each donor based on total available plasma. This allows an overall average increase in plasma volume collected per donation. Unlike the simplified nomogram in place in the United States since 1992, which was based solely on weight, the Persona Nomogram is based on body mass index (BMI) and hematocrit.

#### Persona Nomogram

WEIGHT

HEIGHT


BMI

HEMATOCRIT

**23K**

THE NEXSYS PCS SYSTEM WITH PERSONA TECHNOLOGY, USING THE PERSONA NOMOGRAM, WAS SHOWN TO BE SAFE AND EFFECTIVE IN A CLINICAL TRIAL INVOLVING 23,137 PLASMA DONATIONS—ONE OF THE LARGEST EVER CONDUCTED FOR PLASMAPHERESIS.<sup>1</sup>

### THE PERSONA™ NOMOGRAM

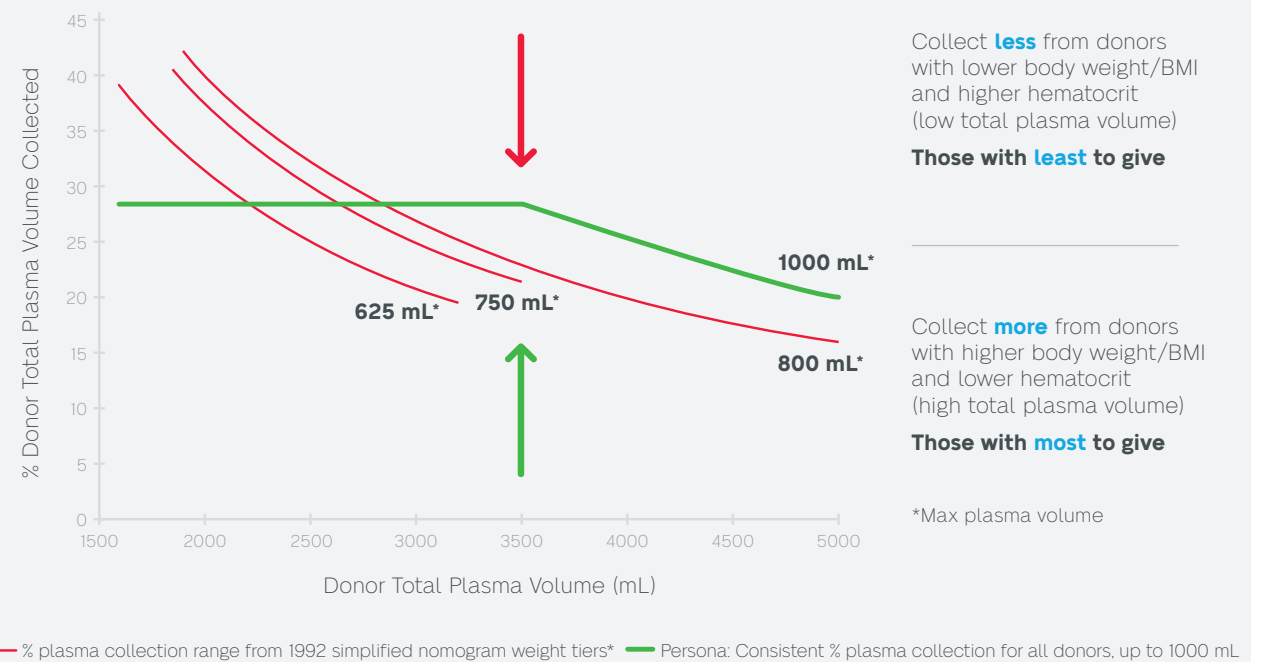

\*FDA Memorandum: Volume Limits – Automated Collection of Source Plasma (11/4/92). Center for Biologics Evaluations and Research. US Food and Drug Administration.



More personal. More plasma.

# Donors come in every size.

Now, so do collection volumes, with the **Persona™ Plasma Collection Solution.**

**Persona™**, a proprietary integrated plasma collection solution built upon the **NexSys PCS® platform**, is the first and only donor-tailored solution clinically shown to yield +9% to 12%\* more plasma per donation on average to maximize both cost-efficient output and patient impact.



**+9% to 12%\*** AVERAGE INCREASE IN PLASMA VOLUME

\*Based on baseline device, software configuration and donor population

### INDIVIDUAL DONOR-TAILORED COLLECTIONS

A **transformational approach** to plasma collection, the Persona™ Plasma Collection Solution enables a more tailored collection target using donors' individual characteristics.

Plasma collection becomes **more personalized** to better align with an individual's body composition.



**NexLynk DMS® software with Persona™ Technology** seamlessly enables plasma collections using the Persona Nomogram.

### HIGHER PLASMA YIELDS

Uses the Persona™ Nomogram to collect a consistent percentage of plasma volume based on height, weight, BMI and hematocrit across all donors, up to a maximum of 1,000 mL, resulting in **an average increase in plasma volume collected per donation†**



**New larger Persona Plasma Pooling Bottle** accommodates increased donation volumes.

### IMMEDIATE AND LONG-TERM GROWTH VALUE

Persona™ can **immediately increase collection volume output** from existing plasma centers, providing a new option to meet the growing demand for plasma. Persona can also increase productivity of potential future plasma center expansions.

Longer-term, more personalized experiences can meaningfully increase loyalty and trust, according to research.†† This potentially leads to **improved quality and longevity of donor relationships**, meeting head-on the difficulty of cost-effectively expanding a donor base and improving retention.



Persona Service and Support offers **fully supported implementation** and **overnight conversion** to minimize disruption to center operations.


†Demonstrated by the IMPACT study, a multicenter, prospective, double-blinded, randomized controlled clinical trial involving 23,137 plasma donations

††Coelho, P. S., & Henseler, J. Creating customer loyalty through service customization. *European Journal of Marketing*. 2012;46(3-4):331-356.

### A TRANSFORMATIONAL APPROACH LEVERAGING THE PROVEN NEXSYS PCS PLATFORM

The Persona™ Plasma Collection Solution is built upon the NexSys PCS®/NexLynk DMS® integrated platform, which has been shown to deliver:

- ↑ plasma yield: +3% (with YES® Yield Enhancing Solution)
- ↑ productivity: increased average turns/device/day
- ↑ quality/compliance: fewer errors
- ↑ donor satisfaction: quieter with more procedure progress information

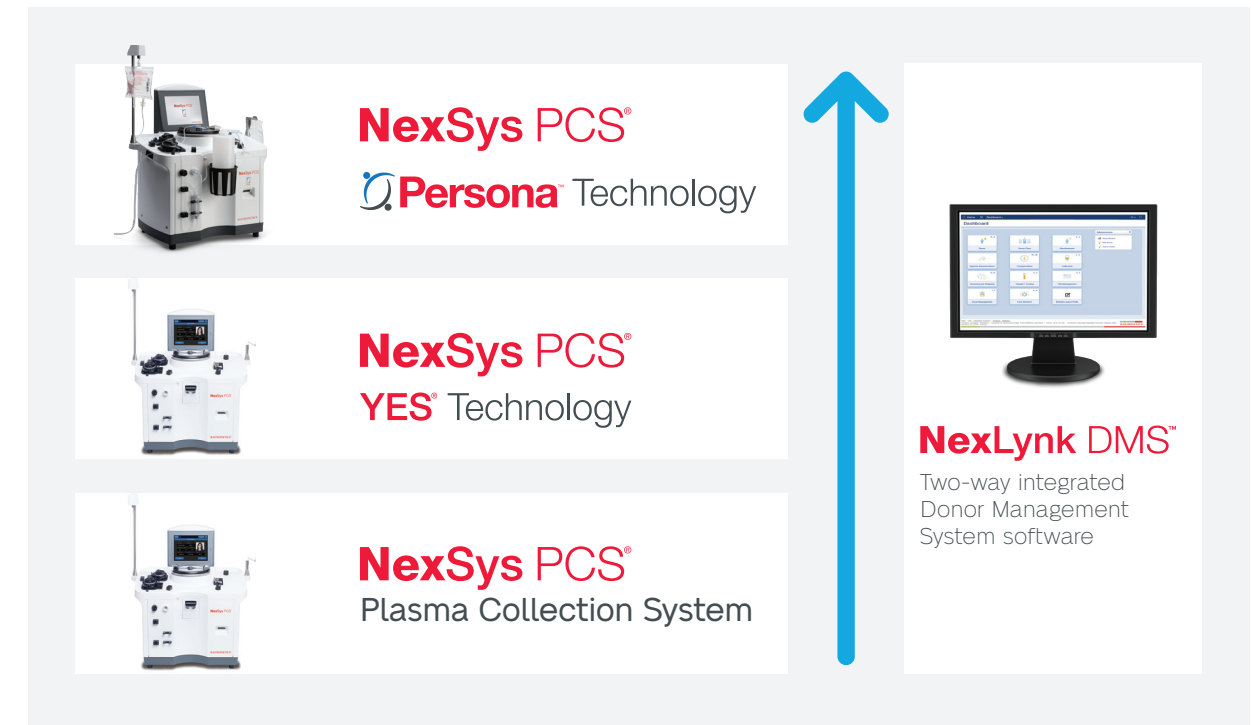

**The world needs more plasma. You can deliver more, with Persona.**

**Contact us today to learn more about how Haemonetics can help tailor your plasma donations and improve your plasma collection productivity.**

For more information on how the novel, proprietary and personalized Persona™ Plasma Collection Solution can create immediate and long-term value for you, contact your local Haemonetics representative today or visit [persona.haemonetics.com](https://www.persona.haemonetics.com).

#### References

1. (2020), Plenary Oral Abstract Session. Transfusion, 60: 7A–260A. <https://doi.org/10.1111/trf.16041>

**For a list of worldwide office locations and contact information, visit [www.haemonetics.com/officelocations](https://www.haemonetics.com/officelocations)**

U.S. Patent Numbers 10,792,416, 10,758,652

© 2020 Haemonetics Corporation. Haemonetics, NexLynk DMS, NexSys PCS, Persona and YES are trademarks or registered trademarks of Haemonetics Corporation in the USA, other countries, or both. 10.2020 USA. COL-COPY-001457-US(AB)

**Terumo v. Haemonetics - PGR2025-00078**  
**Petitioner Terumo Exhibit 1016-0004**